


Product name	Charge massage rod
Product model	CC-9059
Rated voltage	8.4V
Rated power	24W
Rated time	1.5-2 hours
Product materials	ABS PP PA66
Standard of execution	GB4706.1-2005 GB4706.10-2008

Technical parameter



Use operation

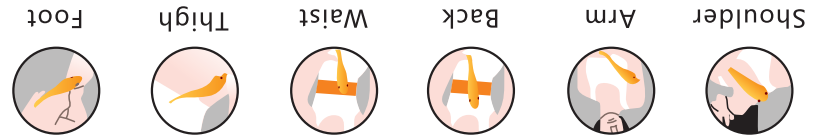

Purpose



CC-9059

WIRELESS RECHARGEABLE MULTI-FUNCTION MASSAGER



Product instructions

Before using the product, please read all the instructions carefully,
Please also kindly store/keep this user manual for future references.

- 1 Before cleaning, switch off the power supply and pull off the power plug. The daily maintenance is clean, please use the neutral detergent and the cloth to wipe, and the corrosive cleaning agent is not available.
- 2 Although the machine is strong, but still please avoid collision or heavy pressure to avoid hard articles to stab the machine surface, so as not to cause unnecessary failure.
- 3 Please use neutral cleaning agent to wash the machine. Do not use gasoline, oil refining, strong acid or strong alkaline cleaning agent.
- 4 Note: the power supply must be removed when cleaning.

Cleaning and maintenance

- ⊗ Postoperative surgery
- ⊗ The bone marrow is fragile or damaged, especially the bone marrow of the back.
- ⊗ Suffering from strain back disease
- ⊗ Suffering from heart disease
- ⊗ Implantation of a cardiac pacemaker or similar implantable medical equipment

In the case of the following, consult your doctor before using this product:

Warning

When using this product, sometimes it will disturb the TV and so on. Please keep the appropriate distance.

Reminder:

1. When the electricity is full, the charger is turned from red to green, that is, it is full of electricity, and the charger is pulled off in time.
2. When the charging plug is inserted into the charging port, the machine can not operate, and the plug can be used.
3. Please do not use in the bathroom, such as high humidity, do not immerse into the water or near the water, to avoid short circuit or electric shock.
4. The power line must be removed after the cleaning and use of the fault and after the end of the use.
5. In lightning, please do not connect the power supply, so as not to cause excessive current, burn down the machine, and even lead to adverse consequences.
6. Verify that the plug is inserted and then used again, otherwise the internal start of the product may be bad.
8. Don't pull the power line directly with your hand when you plug in the power plug.
9. Plug and pull the plug for at least 5 seconds, do not frequently insert or pull.
10. If the power cord is damaged, it must be replaced with a professional soft wire or a special component purchased from its manufacturer or maintenance department.
11. If the internal parts are in trouble, please send them to the service stations all over the company. The company will provide after-sales service for hot pillow.
12. Children and mentally retarded persons with disabilities should be used under the supervision of people.
13. If abnormal noise, odor or temperature are abnormal during work, the use should be stopped immediately.
14. When the machine is not used, be placed in a place that is not easily accessible to children and place children to play or operate machines.

Basic safety specification