



WHOLE NEW EXPERIENCE
360 Degrees All Directional Foot Massaging



CENOCCO CC-9080
The Next Generation Foot Massager
For Blessing, For Perfection, For Health

***Dredging Meridians, Promoting Qi and Circulating Blood Clearing Toxins,
Caring Skin and Relieving Fatigue***

Thank you for purchasing our foot massager!

Please read the user manual carefully before using and pay special attention to safety instructions for proper operation and use of the product

This provides users with safe and efficient approaches of operation, so please properly retain it for future reference and inquiries after reading.

Table of Contents:

- Safety Instructions
- Warranty
- Product Components
- Feature and Use Diagram
- Settings
- Specifications
- Troubleshooting
- Maintenance

Maintenance:

A.) Storage

- The product should be put into the color box and stored under safe, dry, well-ventilated, and cool condition when it is not use.
- Do not twist switching the power button and the cord.
- Do not hang the product with the cord.
- Avoids scratching by sharp objects on the surface of the product.
- Do not put the product in direct sunlight or under condition of high temperature.

B.) Cleaning

- Make sure the power button is off and unplug the power cord before cleaning.
- Use a dry cloth to clean the surface of the product.
- Do not hang the product with the cord

Note: Do not use corrosive detergents, such as benzene or thinner, to clean the product.

C.) Maintenance

- Keep the product clean and away from high humidity conditions.
- Please use a dust cover (not included in the box), to cover the product to prevent dust from coming into it if it is not used for a long period of time.

Troubleshooting:

Please refer to the following several types and possible solutions if problems occur during usage of this product. The product may be returned to the service center or undergone inspection and repair by authorized technicians.

Problems	Possible Causes	Remedies
Not working	Unplugged	Connect the product to power or check whether the product is well-connected
	Power button not pushed	Push the power button
Suddenly stops running	Power disconnection	Check whether the product is well-connected
	Operation duration over 15 minutes	Restart the product
	Product overrun and overheat	Restart the product after cooling
	Too much operating pressure that the massage balls stop rotation	Restart the product

Safety Instructions:

A.) Important Safety Instructions

- Be sure to read the user manual in a detailed way prior to use to ensure safety and proper use of the product
- Make sure to use the product by following the instructions in tis manual
- After using or before cleaning the product, please disconnect the power.
- Always cautiously use the heating function to avoid serious burns.
- Do not apply the product to too insensitive and over sensitive skins or to body parts of poor blood circulations.
- Never allow children or those of physical and mental disorders to use the heating function or improper operation may be dangerous.
- To prevent fire, electric shock or injuries caused by overheating, never use it under a pillow or a blanket.
- Do not cover it during using to avoid malfunction or damage resulting from motor overheating.
- Children, those of physical and mental disorder and the disabled should be under guidance when using the product
- Please stop using immediately if you feel uncomfortable.
- When the power cords or plugs are damaged or dropped into water, do not use this product. Please contact the service center for inspection and repair,
- Do not use spare parts that are not recommended by the manufacturers.
- Do not insert or drop objects into any opening of the products.
- Please close all operational controls before disconnecting the product to the power if you want to cut off the power.
- This product cannot be used for self-diagnosis or treatment and please visit professional medical care if in need of diagnosis or treatment.

B.) Use Environment

- Do not use this product at room temperature over 40°C
- Do not use the product outdoors.
- Do not use the product next to an electric heater or stoves or in direct sunlight.

Safety Instructions:

C.) People Not Recommended Using the Product

- *Please consult your doctor before using if you have the following situations,*
- *Receiving medical treatment or experiencing physical discomfort*
- *Suffering from hypertension, heart disease, cerebral atherosclerosis, and cerebral thrombosis*
- *Suffering from cancer*
- *Pregnant or in a period*
- *Suffering from osteoporosis*
- *Using a pacemaker or implantable medical devices*
- *Suffering from acute diseases*
- *Minors*







D.) Buttons

- *Do not operate the buttons with wet hands.*
- *Keep the buttons from water or other liquid.*
- *Do not fall asleep when using the product.*
- *Do not use this product when you are drunk or feeling uncomfortable.*
- *Do not directly contact the massager pad with your hands, feet joints, knees, chest, and abdomen.*
- *Do not use this product over thirty minutes in a row.*
- *Do not use if the product's one hour before and after eating to avoid adverse effects.*
- *Do not enable the heating function when the product is being used in a low room temperature, it is recommended that you should gradually adjust the temperature to a normal range before enabling heating.*
- *If the product is used immediately after being transferred from a cold environment to a warm environment, the accumulation of water vapor inside the product may affect the normal operation of the product even cause mechanical failure. You are suggested to put the products in a normal room temperature for 1 hour before starting to use it after the transfer*

Specifications:

Name	Airbag Heating Foot Massager
Rated Voltage	220V
Rated Power	50W
Standards	GB 4706.1-2005
Standards	GB 4706.1-2008
Timer	15 minutes

Displays:

	Display
	Timer Button: To set massage duration on timer.
	Mode Button: To switch between massage modes. This product provides two modes of massage with red and blue light indicators for each.
	Intensity Button: To control the massage intensity, where the red light indicates a high intensity and blue light indicates a low intensity
	Heating Button: When the product is turned on, push this button to control the heating independently with a red-light indicator (The heating is on is default)
	Power Button: To turn on/off the product

A.) Before use

- Check the product for any damage.
- Ensure a stable and comfortable use environment.

B.) Use

- Untie the power cord and connect to power of 220V/50Hz


Default Settings

- Press the Power Button Once, and the heating function is ON in a red indicator: The Intensity is low in a blue indicator, and the product is on Mode 1 in a blue indicator.
- Press the Power Button ONCE, then press the Heating Button ONCE, and the heating function will be OFF as the red indicator goes off.
- When the product is in operation, press the Intensity Button to switch between the two types of intensity. Blue light indicates low Intensity and red light indicates high Intensity.
- When the products is in operations, press the MODE Button to switch between two modes. Blue light indicates Mode 1 and red light indicates Mode 2.

Warranty:

1. Warranty is not transferable.
2. Users can be enjoying free warranty from the date of purchase during the warranty period) The surface is not covered under warranty.)
3. Product warranty is valid for one year.
4. Please present the original proof of purchase and warranty card when applying for warranty service.
5. The purchasing date shall be based on the effective receipt or invoice issued by the retailers.
6. This warranty shall not apply to leasing, trade, commence, institutions and other non-buyer users.
7. The following situations are not covered by the warranty and only paid service are provided (charging labor costs and materials costs)
 - Destruction caused by human factors, including the non-normal using environment, usage and storage that are not complaint to the instructions, and damage caused by improper use. etc.
 - Disassemble, repair and modification of the products without our authorization.
 - Damage caused by user poor transportation after purchase.
 - Damage resulting from any other force majeure (such as flood, fire, lighting, and earthquakes).

Product Components:

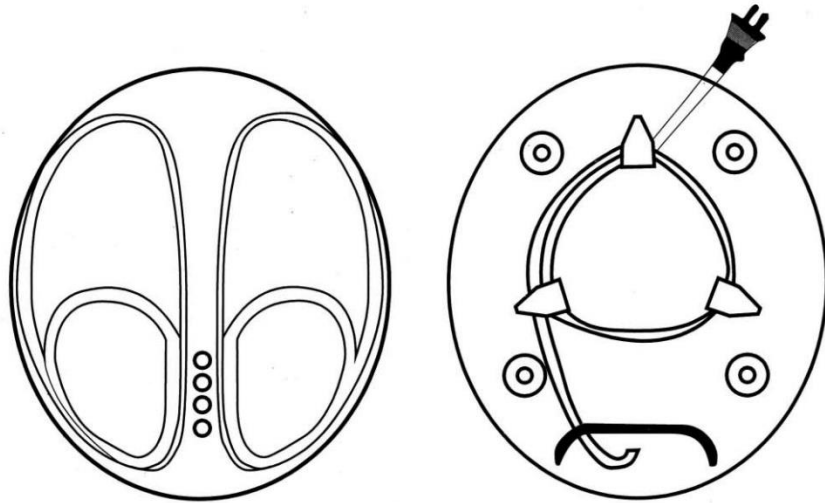

1. **Central Surface: Four Independent Buttons** - Controlling the product's settings and functions
2. **Central Massager Pad: Scratching** - Acupressure on the arches with heating.
3. **Front Massager Pad: Roller Kneading** - Acupressure on the front soles.
4. **Rear Massager Pad: Push Kneading** - Acupressure on the heels
5. **Interior airbag** - Air pressure on the insteps
6. **Removable and washable fabric cover** - cleaning
7. **Bottom: Cord Wrap** - keeps things tidy

Features and Use Diagram:

Product Features:

- *Egg-shaped design, fashionable and beautiful, lightweight, and portable:*
- *A combination of a variety of functions including scratching, roller kneading and push kneading acupressure, air pressure and heating.*
- *Massager Pad that provides all-around massages and in-deep acupressure to clear the body's meridians and foot reactions zones.*
- *Removable and washable cover, allowing you to keep clean at any moment.*
- *Cord Wrap, bringing convenience to life.*

Use Diagram:

