



EMS-EYE MASSAGER CC-9099

User Manual

Thank you for buying this product.
Read the instructions carefully before using the product.
Please keep the instructions properly.

PRODUCT INTRODUCTION

EMS-Eye Massager is a high-end eye beauty instrument with EMS Micro Current+ Beauty Red Light+ Constant Temperature Hot Compress+ High Frequency Vibration and other functions. Long term use can effectively remove the black eye bags, relieve the eye skin problems such as sagging, corner lines, dark spots, etc... effectively restore skin elasticity and delay skin aging.

INSTRUCTIONS

1. Press and hold the key to turn on/off. and the default Its mode is release mode.
2. Double click the key to select the mode according to your needs.

A. Relaxation Mode: High frequency vibration+ 42°C Constant temperature hot compress, promote blood circulation of eyes. relieve fatigue, remove eye bags, and dark circles; accelerate the absorption of eye cream and other skin care products.

B. Beauty Mode: High frequency vibration+ Beauty red light, 620nm beauty red light for rough and dark skin, can smooth microcirculation, rejuvenate skin energy, and effectively fade freckles and fine lines.

C. EMS Tightening Mode, EMS+ Red light. EMS direct to the dermis through micro current simulation cell bioelectricity to promote collagen production, promote lymph detoxification. tighten the skin. make the skin full and tight. and effectively relieve the eye skin loose and drooping.

3. Click the key to select the intensity gear. There are 3 gears in total.

4. Please charge for about 2 hours to full charge when using for the first time. When charging, the power display icon flashes; when the power is low, the power display icon flashes for 5 times. and then the power is turned off: when the power is full , the power display icon stays on.


EMS

EMS, acts on the muscles of subcutaneous tissue. By simulating cell bioelectricity with micro current. facial muscles are cut to produce more ATP (ATP is the energy needed for collagen synthesis). To increase collagen needed for collagen synthesis). to increase collagen. tighten the skin and promote lymphatic detoxification.

BEAUTY RED LIGHT

Beauty red light is called bioactive light. It transforms light energy into energy in eel ls. accelerates eel l growth process circulation. stimulates fiber cells to produce collagen.so as to increase skin elasticity. improve skin yellowing and dullness and bring unexpected effects on repairing aging skin. acne skin. spot skin. wrinkled skin. allergic symptoms, sunburn skin. And will not cause any harm and discomfort to the skin.

FUNCTIONS

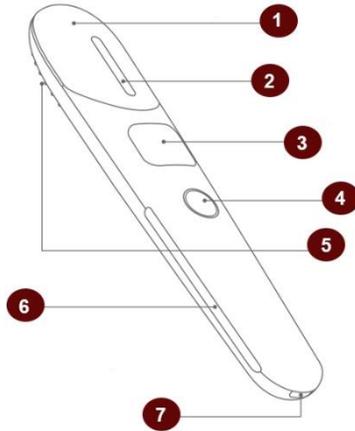

(1) Metal Working Head	(2) Beauty Red Light
(3) Display	(4) Key
(5) Massage Dot	(6) EMS Conduction Decoration
(7) Charger Port	

EMS CONDUCTION DECORATION

(It is necessary for the skin of the hand to contact the conductive strip, so that the metal working head can generate EMS micro current)

WORKING INSTRUCTIONS

Double click the key conversion mode, and click the key to convert the strong master gear:

♥ RELAXATION MODE	(O) + Hot
♥ BEAUTY MODE	(O) + Light
♥ EMS COMPACT MODE	EMS + Light

EMS intensity is divided into 3 levels. Please select the appropriate intensity according to your skin sensitivity to avoid discomfort. It is recommended to start from low gear:

	1 light display. EMS intensity is “Weak”
	2 light displays, EMS intensity is “Medium”
	3 light display, EMS intensity is “Strong”

SPECIFICATIONS

Name	EMS-Eye Massager
Red Band	620±20nm
Heating Temperature	42°C±2°C
EMS Current	Up to 30µA
Rate Voltage	3.7V
Heating Power	1.9W
Rate Current	450mA
Battery Capacity	500mAh
Work Time	More than 90 minutes
Charging Time	About 2 hours
Product Size	160 x 30 x 19.5mm
Executive Standard	GB4706.1-2005 GB4706.15-2008

CLEANING METHOD

After the product is used, wipe it off with a cloth in the off state and put it in a cool and dry place. Note that this product is not waterproof and should not be rinsed with water.

ATTENTION

1. This product has no waterproof function. please do not put it in the water.
2. The product has heating function, about 42°C. People who are not sensitive to heat should use it carefully.
3. When the eyes are painted foundation or make-up, please remove the makeup before taking care.
4. If you feel discomfort or skin abnormality during use, please stop using immediately and consult your doctor.
5. People who are not allowed to use this product:

A. The corresponding parts with metal in the body due to treatment and other reasons:

B. Patients with heart disease (especially those with pacemakers)

C. Capillary expansion caused by hormone or disease.

D. Skin inflammation, pruritus, redness and swelling caused by allergy or other reasons.

E. It is not recommended to use this product for the parts with micro finishing, prosthesis, nose augmentation, injection, embedded gold thread, etc.

6. Do not let children touch or use this product to avoid danger.



